

PC Cloning Software for Portable Particle Counters

Backup/Restore, Management & Control



When Accuracy Matters!™

Easy control and management.

Backup & restore of configuration settings.

Intuitive graphical user interface (GUI).

Fast and reliable standardization of portable particle counters.

Computer Requirements

1. Windows XP, Vista, Win 7, Win 8
2. Only requires 10 MB disc space
3. Any Windows compatible computer with either an RS-232 port or USB port

A unique and standard feature of all Climet portable particle counters is the ability to clone settings between two units. This **Unit-to-Unit Cloning** copies programs, location IDs, user IDs, alarm settings, report settings, security settings, and other setup parameters from one particle counter and copies to another.

For users with more complex business requirements, Climet offers our **PC Cloning Utility**. This software solution permits users to upload a Climet particle counter's configuration to a PC; or alternatively, users can create custom settings on their computer that can be downloaded to their Climet particle counters.

The advantage of PC Cloning Utility is the ease of use when creating, editing or copying configuration settings. All settings and programs can be easily managed from a central PC and then seamlessly uploaded to a Climet Particle Counter with the simple click of the mouse. This provides a trouble-free means of ensuring multiple particle counters have a standard configuration that meets your security and monitoring requirements.

Applications

The PC Cloning Utility is a perfect solution for companies with multiple particle counters that need to be configured identically, or managed centrally. Additionally, companies with multiple users of a single particle counter will use PC Cloning software to easily manage users and administrative privileges. Cleanroom Certifiers also use PC Cloning when custom IDs are needed for multiple customers. Finally, it is a perfect solution for easy backup and restore of particle counter configuration settings.

The software utility comes complete with a CD-ROM and an RS-232 serial adapter with USB cable for connecting the particle counter to either another particle counter, or to your laptop, workstation, or server.

Applications Include:




Cleanroom Monitoring & Validation
Medical Device Manufacturing
Pharmaceutical Manufacturing
Hospitals
Semiconductor & Disk Drive Clean

Food & Beverage Processing
Aerospace Assembly
room Monitoring
Cosmetic Manufacturing
Research & Development

CLIMET[®]
INSTRUMENTS COMPANY

PC Cloning Software

Part Number	95135100
Compatibility	Compatible with all Climet CI-x5x series of portable particle counters with v 3.00 firmware or later.
Standardization	Ideal solution when two or more portable particle counters are used that need identical security and other monitoring parameters.
Backup & Restore	Conveniently stores Climet portable particle counter settings and configuration parameters. Ensures a trouble-free restoration of these settings if they are ever lost or deleted for any reason.
Special Applications	<p>Ideal when multiple users are operating a single particle counter. Greatly simplifies administration of User IDs and other configuration parameters.</p> <p>Used when more than 30 programs or more than 100 location IDs are required.</p> <p>Used for configuration of instruments for certification, then convert back to monitoring.</p>
Included Accessories	<p>CD-ROM</p> <p>DB-9(f) to DB-9(f) cable.</p> <p>DB-9(m) to USB Multiplexor Adapter Cable. (converts RS-232 to USB protocols)</p>
Hardware Requirements	<ol style="list-style-type: none"> 1. Window XP, Vista, Windows 7, and Windows 8 2. Windows compatible laptop, workstation, or server 3. 10 MB Disc Space

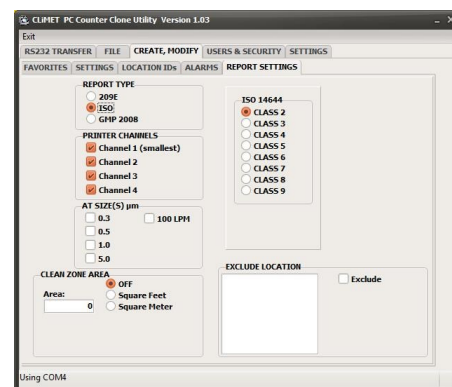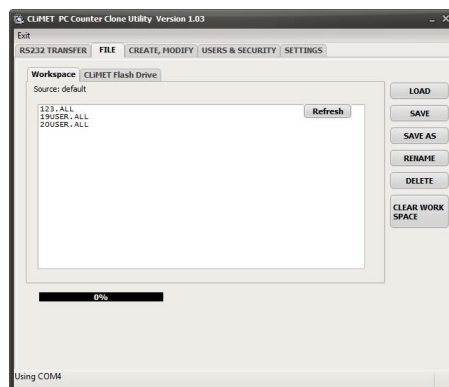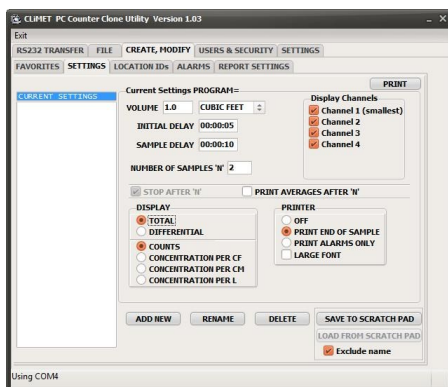

For further technical information and questions, please contact us. **We value the opportunity to work with you.**

Climet Instruments Company

1320 West Colton Ave
 Redlands, CA 92374 (USA)
 Tel: +1 (909) 793-2788

www.climet.com | sales@climet.com

