



CLIMET INSTRUMENTS COMPANY

DataPro 2.5e Software For Portable Particle Counters

21 CFR 11 Compliant

When Accuracy Matters!™



Intuitive Graphical User Interface

Sample Data Viewing

Sample Data Management

Scheduling

Remote Control

Graphing & Trend Analysis

Reporting ISO 14644-1, EU GMP Annex 1, and FS 209E

Compliant with 21 CFR Part 11

Auditing Features

Backup/Restore supported in Windows operating system.

Computer Requirements

1. Windows XP, Vista, Win 7, Win 8, Win 10
2. RS-232 COM Port, or USB to RS-232 converter (p/n: 77308100) needed.
3. Any Windows compatible computer

DataPro 2.5e Software is used with Climet portable particle counters. It provides an intuitive graphical user interface (GUI) for sample data viewing, data management, scheduling, remote control, graphing and trend analysis, as well as reporting and auditing.



DataPro 2.5e implements a robust and comprehensive security protocol, including sample data authentication that is compliant with **FDA 21 CFR Part 11**. It enables our customers to sustain current industry “best manufacturing practices” (CGMP). The software’s integrated audit feature eases compliance and monitoring analysis, provides for faster validation, and better CGMP process management.

This comprehensive software suite provides pass/fail conformance reports to ISO 14644-1, EU GMP Annex 1, and the FS 209E standards.

Graphing and trend analysis is a snap. With the simple click of the mouse, sample data is automatically and seamlessly imported into a graph in Microsoft Excel without having to copy or paste data. DataPro 2.5e provides users the ability to easily generate graphs and other analytics key for trend analysis and early failure identification.

The software suite also comes with Climet’s IDPro utility at no additional charge, which is used for creating and storing IDs on a workstation, server or laptop.

Applications

Designed for pharmaceutical, biotechnology and other companies that do not have a Laboratory Information Management System (LIMS). This is a perfect solution for those organizations seeking to deploy an affordable high-reliability secure facility monitoring suite for their Climet portable particle counters.

Applications Include:


Cleanroom Monitoring & Validation
 Medical Device Manufacturing
 Pharmaceutical Manufacturing
 Hospitals
 Semiconductor & Disk Drive Clean

Food & Beverage Processing
 Aerospace Assembly
 room Monitoring
 Cosmetic Manufacturing
 Research & Development




DataPro 2.5e Software

Part Number	95037000
Compatibility	Compatible with Climet CI-x5x Series and CI-x70 Series of portable particle counters.

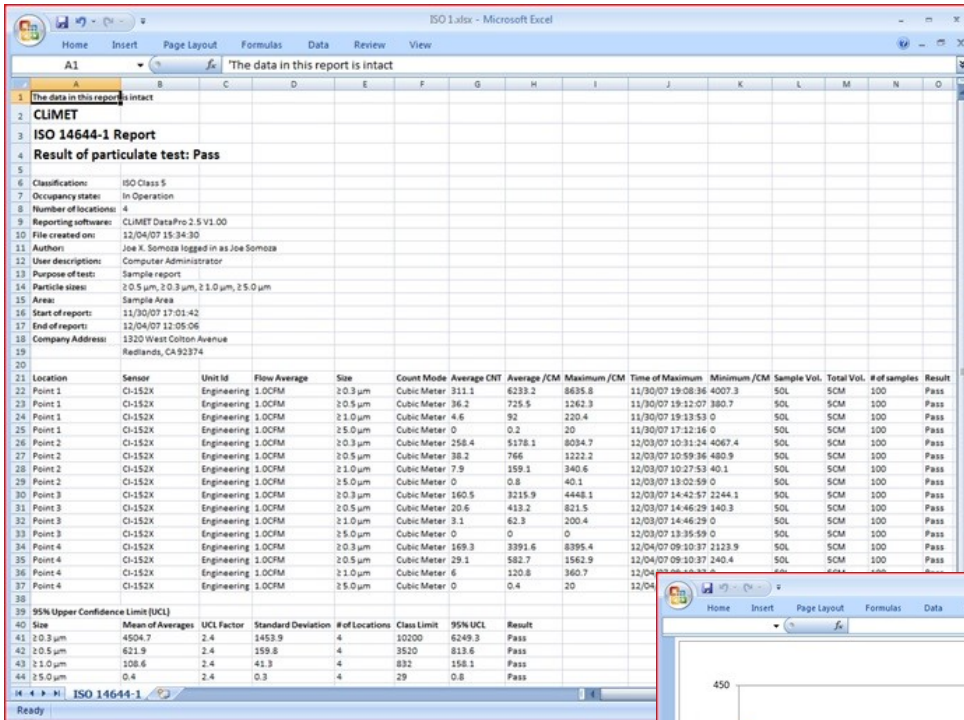

To the left is a sample conformance report for ISO 14644-1. Other reports are available for EU GMP Annex 1, and FS 209E.

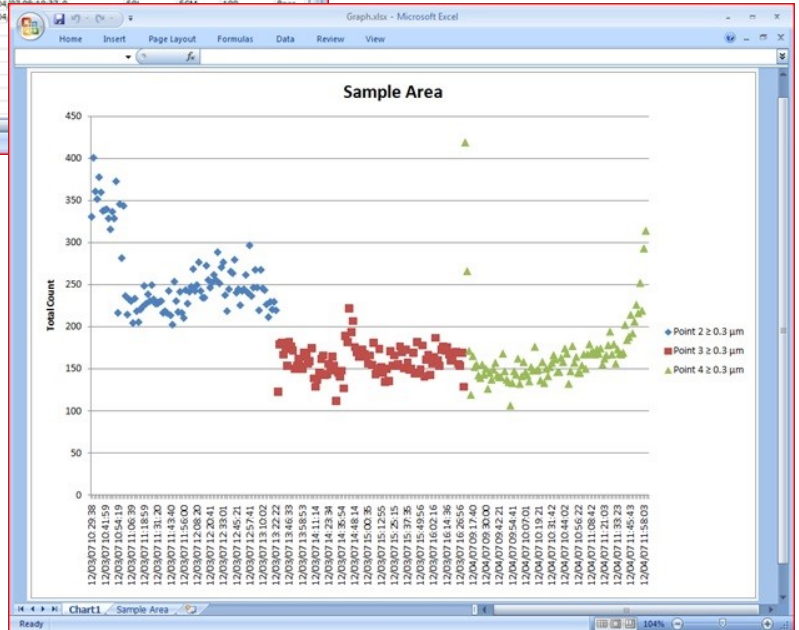

To the right is a graph used for trend analysis.

For further technical information and questions, please contact us. **We value the opportunity to work with you.**

Climet Instruments Company

1320 West Colton Ave
Redlands, CA 92374 (USA)

Tel: +1 (909) 793-2788

www.climet.com | sales@climet.com

